

Otto Bock[®]

QUALITY FOR LIFE



Manual Wheelchairs

Wheelchairs for Active Use | Standard Adaptive Wheelchairs



Wheelchairs for Active Use



Standard Adaptive Wheelchairs



Frame Colours

Table of Contents

Wheelchairs for Active Use	4
Blizzard	6
Voyager	12
Avantgarde T	18
Avantgarde Ti 8.9	24
Avantgarde Ti Ultra	30
Avantgarde VR	36
Avantgarde XXL 2	42

Standard Adaptive Wheelchairs	48
StartM2 Effect	50
StartM3 Hemi	56
StartM4 XXL	62
StartM5 Comfort	68
StartM6 Junior	74

Frame Colours	80
----------------------	-----------

The Blizzard sustains its excellent reputation due to its reliability and suitability for use in all circumstances. It takes only a few easy steps to lay the wheelchair compactly on the passenger's seat. Simultaneously its rigid frame, which is fabricated individually and offers ideal geometry for its user, provides safe pleasure while driving. In short: The Blizzard combines an attractive and award-winning design with functionality and a sporty lifestyle. A real eye-catcher!

The product range in the rigid-frame sector is complemented by the proven Voyager with its brawny, timeless design. In addition to being very adaptable and capable of being loaded into a car by the user, the Voyager shines due to its extremely torsion-resistant structure – and not just in the polished version.

It can meet any challenge, both on and off the asphalt. A powerful package with character.

The success of the Avantgarde wheelchair series simply lies in the fact that: Active wheelchair users just have more fun thanks to the manoeuvrability, variations and lightness of the models. We have continued to search for further improvements and have now arrived at even lighter, more manoeuvrable and more flexible solutions. Another one of the Avantgarde family's greatest strengths is that its very appealing design meets the increasingly common desire not only to be mobile, but also to amplify individuality and to stand out attractively in doing so.





Wheelchairs for Active Use

With its sporty, elegant design, the Blizzard won the "red dot design award" for high design quality. In addition, the Blizzard offers all the benefits of a wheelchair for active use: optimal handling, a compact folding size for easy loading, and a precisely adjustable centre of gravity. With its broad range of options, the Blizzard is mainly intended for active people who demand an attractive and dynamic wheelchair. For example, the Blizzard can be equipped with a wheelbase extension for integration with a hand bike.



Blizzard





Designer caster with angled caster fork,
for optimal driving behaviour



View into the open frame with rigid axle



Clothing protector side panel with
eccentric disc for fine adjustment



Seat unit with adaptable upholstery
and titanium stabilizer



Folding push handles – available, if needed;
out of the way, if the user him- or herself
propels the wheelchair



Rear wheel with asymmetric
Downtown tyre and titanium
push ring

Blizzard



Options

Footrest tube with footplate

- Anti-slip
- Optional lateral heel blocks
- Angle-adjustable



Adjustable seat upholstery

- Individual adaptation with variable hook and loop system



Folding, lockable back

- Easy transport and compact folding size



Blizzard

Product Features

- Folding back
- Variable seat height
- Variable seat depth
- Variable wheelbase
- Very manoeuvrable
- Rigid-frame construction



Technical Data

Seat width	32 cm - 46 cm	12.6 in - 18.1 in
Seat depth	36 cm - 48 cm	14.2 in - 18.9 in
Back height	25 cm - 50 cm	9.8 in - 19.7 in
Weight	from 9.5 kg	from 20.9 lbs
Front seat height	44 cm - 53 cm	17.3 in - 20.9 in
Rear seat height	38.5 cm - 54.5 cm	15.2 in - 21.5 in
Lower leg length	17 cm - 52 cm	6.7 in - 20.5 in
Max. user weight	125 kg	275 lbs
Front frame angle	75° / 85°	75° / 85°

Robust and indestructible – with the Voyager, the name is the game. During active everyday life and when traveling, this agile all-rounder excels with precise handling and extremely good stability. The compact design with an open frame and folding backrest makes it possible for the user to load the wheelchair into a car by him or herself. Moreover, the Voyager is easy to manoeuvre in restricted spaces.

Active users often drive with an inclined seat. Therefore, the axle geometry and backrest angle adjustment allow for extreme settings in order to turn the seat and back into a single ergonomic unit. Various options expand the range of applications, while the Voyager Light Edition is available for purists. Reduced to the essentials, the smallest size of the Voyager has a total weight of about 9 kg. Ultra-light wheels, aluminium screws, breathable upholstery, and a titanium tube footrest are standard equipment.



Voyager





Luminous wheels provide fun and passive safety



The open tube footrest offers maximum stability and is included with the Light Edition



The Infinity Ultralight rim noticeably reduces the total weight of the Voyager



The compact wheel lock with short operating lever is not in the way during transfer



Adaptable back upholstery for optimal setting as desired



The clothing protector reliably prevents clothing from getting dirty

Voyager



Options

Shock absorber

- Absorbs shocks caused by uneven surfaces
- Prevents the occurrence of spasms
- Different degrees of stiffness available



Side panel with armrest

- Offers enhanced support for the arms
- Height-adjustable without the use of tools



Footplate with lateral heel blocks

- Offers maximum space for the feet
- Different angle and depth settings are possible
- Lateral heel blocks prevent the feet from slipping off the footplate



Voyager

Product Features

- Open frame construction
- Folding, angle-adjustable backrest
- Continuously adjustable centre of gravity
- Rear seat height adjustable in small increments of 1 cm
- Offset of 0, 2 or 4 cm can be selected
- Continuous adjustment of lower leg length
- Optional shock absorber system



Technical Data

Seat width	34 cm - 46 cm	13.4 in - 18.1 in
Seat depth	36 cm - 48 cm	14.2 in - 18.9 in
Back height	25 cm - 45 cm	9.8 in - 17.7 in
Weight	from approx. 9.8 kg	from 21.6 lbs
Front seat height	43.3 cm - 53 cm	17.1 in - 20.9 in
Rear seat height	40 cm - 55.5 cm	15.7 in - 21.8 in
Width	66.5 cm - 82 cm	26.2 in - 32.3 in
Back angle	-15° - +27°	-15° - +27°
Front frame angle	75°, 85°	75°, 85°
Lower leg length	17 cm - 52 cm	6.7 in - 20.5 in
Max. user weight	120 kg	265 lbs

The Avantgarde T is the wheelchair for all occasions. It is individually configured from an extensive assembly kit and represents the basis of the entire Avantgarde family. Functionality has number-one priority in the standard model. The Avantgarde T is especially stable, thanks to its firm aluminium frame with a double crossbrace. Moreover, the Avantgarde T features outstanding manoeuvrability and a vast range of options such as the optional Ergo back.



Avantgarde T






Business upholstery with various practical storage compartments for purse, mobile phone, palmtop, etc.



Push handle with non-slip grip




Single-panel footrest with angle-adjustable aluminium footplate

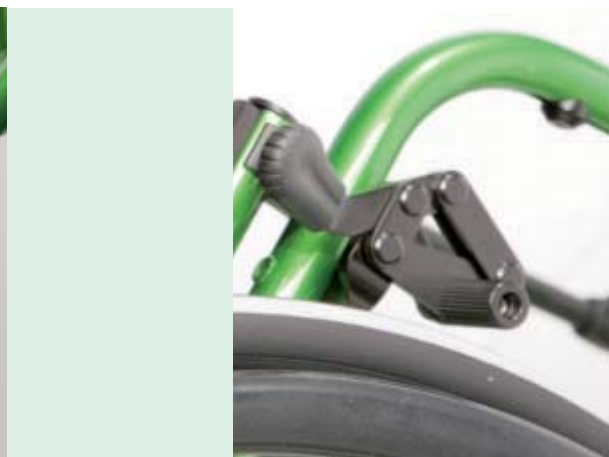

The wheel lock with short operating lever is not in the way during transfer



The plug-on clothing protector is easy to remove

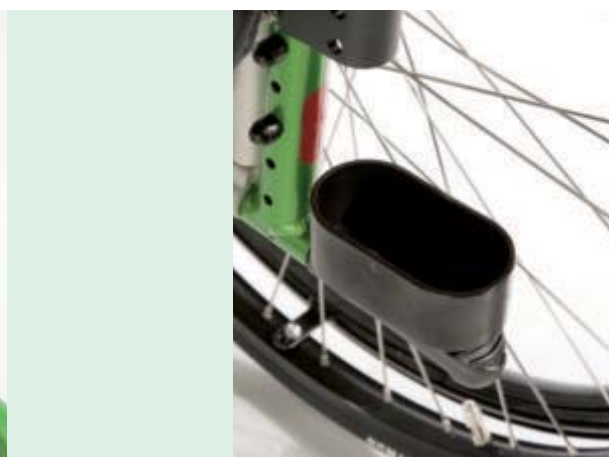

The crutch holder offers space for two forearm crutches

Avantgarde T



Options


Removable tube armrest

- Padded
- Very light



Shock absorber

- Rear wheel suspension for high driving comfort

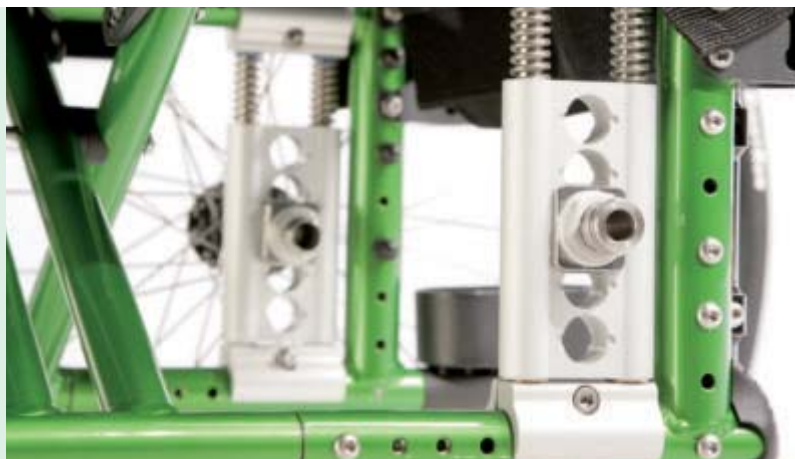

Push ring cover

- Better grip
- Enhanced force transmission



Avantgarde T

Product Features

- High-strength aluminium frame with double cross-brace
- Single-panel footrest or individual footrests with "easy click" locking function
- Height-adjustable and adaptable back as standard equipment
- Aluminium side panels
- Precise setting of the rear wheel position by means of optimized adapter
- Front frame offset can be selected
- A vast range of options such as Ergo back with physiologically adjustable pivot point and angle adjustment from -9° to +9°



Technical Data

Seat width	32 cm - 50 cm	12.6 in - 19.7 in
Seat depth	36 cm - 48 cm	14.2 in - 18.9 in
Back height	25 cm - 50 cm	9.8 in - 19.7 in
Front seat height	39 cm - 56 cm	15.3 in - 22 in
Rear seat height	37 cm - 53 cm	14.6 in - 20.9 in
Seat inclination	adjustable up to 15°	adjustable up to 15°
Weight	from approx. 11.5 kg	from approx. 25.3 lbs
Max. user weight	125 kg	275 lbs

The Avantgarde Ti 8.9 attracts attention and is therefore the right model for demanding individualists. It is available as a polished special fabrication or with motifs from the airbrush design line – the creative possibilities are virtually endless with this wheelchair for active use. The Avantgarde Ti 8.9 is especially lightweight and deliberately dispenses with heavy options. Its high-quality ultralight wheels reduce additional weight and increase comfort. Furthermore, the rear wheel position and the centre of gravity of the Ti 8.9 can be adjusted to meet the user's needs. It also features an anatomically adjustable back as well as an angle- and height-adjustable footrest.



Avantgarde Ti 8.9





Adaptable back upholstery, adjustable without the use of tools



Front frame with model description; 8.9 kg is the total weight of the wheelchair



Standard equipment includes side panels with protection against cold made of neoprene



An optional zippered pocket in the back upholstery provides additional storage space

Avantgarde Ti 8.9



Options

Caster mounting with quick-release axle

- Threaded axle is also possible
- Different caster trailing settings are possible



Frame pad

- To protect the legs as well as furniture



Performance wheel lock

- Space-saving and effective



Avantgarde Ti 8.9

Product Features

- Infinity ultra-lightweight rear wheels
- Precise setting of the rear wheel position by means of optimized adapter
- Anatomically adaptable back with strap system
- Breathable back cover
- Aluminium side panels
- Height- and angle-adjustable footrest of aluminium/ titanium components
- Front frame offset can be selected
- Aluminium screws
- Shock absorber system



Technical Data

Seat width	32 cm - 50 cm	12.6 in - 19.7 in
Seat depth	36 cm - 48 cm	14.2 in - 18.9 in
Back height	25 cm - 50 cm	9.8 in - 19.7 in
Front seat height	39 cm - 54.5 cm	15.3 in - 21.5 in
Rear seat height	40 cm - 53 cm	15.7 in - 20.9 in
Seat inclination	adjustable up to 15°	adjustable up to 15°
Weight	from approx. 8.9 kg	from approx. 19.6 lbs
Max. user weight	125 kg	275 lbs

The Avantgarde Ti Ultra is the lightweight among wheelchairs for active use. It is the ideal companion for experienced wheelchair users with an unchanging sitting position. Thanks to its firmly welded rear frame, the Avantgarde Ti Ultra combines two attributes: It offers highest stability and is yet comfortably to transport without requiring a lot of space. A single, firmly welded tube with only one bore hole replaces the obligatory rear wheel adapter. The anterior seat height and the wheelbase can be exactly adjusted to the user's individual requirements. In addition, the modular design of the Ti Ultra allows for combination with many options of the Avantgarde series.



Avantgarde Ti Ultra





The flip-up tube footrest combines function with lightweight construction



Caster fork with caster trailing setting as a harmonious unit



Firmly welded, short push handles



Breathable, adaptable back upholstery



Ultralight aluminium frame; the complete wheelchair is available with a weight starting from 8 kg

Avantgarde Ti Ultra



Options

Side panels

Various versions are possible:

- Carbon
- Long armrest
- Padded
- With clothing protector



Folding push handles

- Ideal for active wheelchair users driving themselves
- With the Ti Ultra available upon request



Closed footplate

- Comfort lock for more stability
- Folds away for easy transfer



Avantgarde Ti Ultra

Product Features

- Infinity ultra-lightweight rear wheels
- Firmly welded rear wheel position
- Anatomically adaptable back with strap system
- Breathable back cover
- Aluminium side panels
- Height- and angle-adjustable footrest of aluminium/ titanium components
- Front frame offset can be selected
- Aluminium screws



Technical Data

Seat width	32 cm - 50 cm	12.6 in - 19.7 in
Seat depth	36 cm - 48 cm	14.2 in - 18.9 in
Back height	25 cm - 50 cm	9.8 in - 17.7 in
Front seat height	39 cm - 54.5 cm	15.3 in - 21.5 in
Rear seat height	40 cm - 53 cm	15.7 in - 20.9 in
Seat inclination	Seat inclination up to 15° can be ordered; cannot be changed later	Seat inclination up to 15° can be ordered; cannot be changed later
Weight	from approx. 8 kg	from approx. 17.6 lbs
Max. user weight	125 kg	275 lbs

The Avantgarde VR is a master of metamorphosis among wheelchairs for active use. Like the Avantgarde T, it offers the full range of options. The Avantgarde VR combines harmonious function and design and is ideally suited for patients with many different indications, such as multiple sclerosis or rheumatic diseases. Transfer of the wheelchair user is safe and easy thanks to footrests that can be swung inward and outward.



Avantgarde VR





The side panel can be swung completely out of the way for easy transfer



Stable double crossbrace, torsion-resistant and easy to fold



Armrest with PU pad, height-adjustable



Adaptable back upholstery, adjustable without the use of tools

Avantgarde VR



Options

Adjustable legrests

- Height-adjustable
- Angle-adjustable
- Swing-away
- Removable



Swing-away anti-tipper

- Different depth settings possible



Tray

- Adjustable
- Tray edge



Avantgarde VR

Product Features

- High-strength aluminium frame with double cross-brace
- Individual footrests with "easy click" locking function, continuously height- and angle-adjustable
- Anatomically adaptable back with specially developed strap system, different height settings are possible
- Stable aluminium side panels as standard equipment
- A vast range of additional options such as Ergo back with physiologically adjustable pivot point and angle adjustment from -9° to +9° and shock absorber system



Technical Data

Seat width	32 cm - 50 cm	12.6 in - 19.7 in
Seat depth	36 cm - 48 cm	14.2 in - 18.9 in
Back height	25 cm - 50 cm	9.8 in - 19.7 in
Front seat height	39 cm - 56 cm	15.3 in - 22 in
Rear seat height	37 cm - 53 cm	14.6 in - 20.9 in
Seat inclination	adjustable up to 15°	adjustable up to 15°
Weight	from approx. 11.8 kg	from approx. 26 lbs
Max. user weight	125 kg	275 lbs

With the Avantgarde XXL 2, it's the size that counts. This wheelchair for active use was specially designed for bariatric patients as well as particularly tall and stout wheelchair users. The Avantgarde XXL 2 has maximum seat width of 62 centimetres and can hold a maximum of 180 kilograms. As a result, the most stable folding wheelchair of the Avantgarde family can withstand maximum stress and strain and still be easily transported without requiring a lot of space.




Avantgarde XXL 2






Back stabilizer, increases overall stability and can be swung out of the way for folding

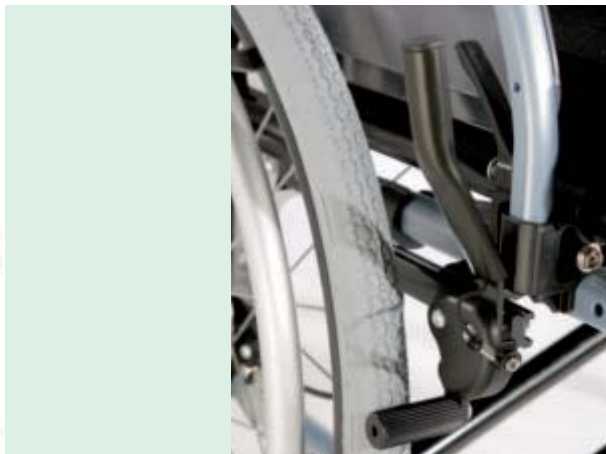

The wheel lock lever extension enables users with limited hand strength to use the wheel lock



Brake lever for attendant

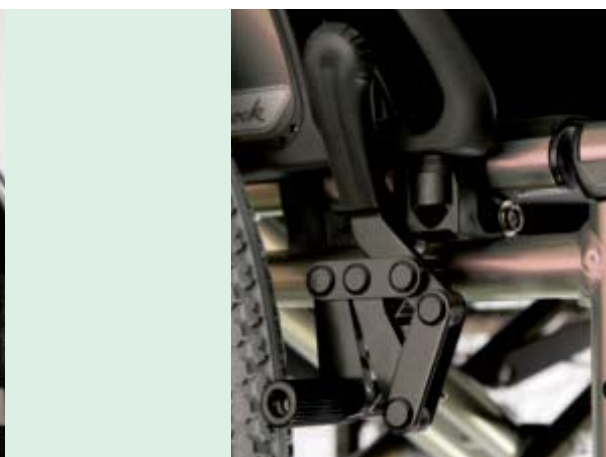

Easy to operate wheel lock in a compact design



Adaptable back upholstery, adjustable without the use of tools

Avantgarde XXL 2



Options

Transport wheels

- Required for small doors



Individual footrests

- Removable
- Angle-adjustable
- Swing-away



Height-adjustable armrest

- Padded
- Long and short armrest



Avantgarde XXL 2

Product Features

- Seat upholstery with reinforced insert
- Light metal frame with vertical front frame part
- Aluminium side panels with neoprene inserts
- Single-panel footrest, angle-adjustable, swing-away and removable
- 24" or 26" rear wheels with quick-release axle, larger hub, hollow rim, and aluminium push ring
- 0° rear wheel adapter
- Low own weight
- Optional back stabilizer



Technical Data

Seat width	46 cm - 62 cm	18.1 in - 24.4 in
Seat depth	40 cm - 56 cm	15.7 in - 22 in
Back height	25 cm - 50 cm	9.8 in - 19.7 in
Variable armrest height	20 cm - 30 cm	7.9 in - 11.8 in
Weight	15.2 kg	33.5 lbs
Front seat height	45 cm - 53.5 cm	17.7 in - 21.1 in
Rear seat height	40 cm - 53.5 cm	15.7 in - 21.1 in
Width	66.5 cm - 82 cm	26.2 in - 32.3 in
Overall height	70 cm - 103.5 cm	27.6 in - 40.7 in
Overall length	79 cm - 114 cm	31.1 in - 44.9 in
Lower leg length	32 cm - 55 cm	12.6 in - 21.6 in
Max. user weight	180 kg	397 lbs

Complex solutions do not have to be complicated. Behind the diversity of the new Start generation is a sophisticated concept whose strengths lie in its straightforwardness and simplicity. Otto Bock, with its new Start series, is now introducing the modular design into the area of standard adaptive wheelchairs.

The result is a high degree of individuality in the fitting, while simultaneously making this easier for specialist dealers. All members of the new Start series are distinguished by their technical-, customer- and user-friendliness, since identical components are used for the entire spectrum of the product family, depending on the target group.

A large variety of product options usually represents a great challenge when it comes to assembly or service. This is not the case with the Start series. The modular design makes it possible for almost exactly the same spare parts to be used in all Start wheelchairs. This significantly reduces stockkeeping. And it also lowers the cost of keeping tools on hand.





Standard Adaptive Wheelchairs

The Start M2 Effect is the purist of the Start product family. This basic model has the same frame module as the more complex versions, but we deliberately chose not to provide excessive options. With its profiled aluminium frame and a uniform caster fork for all wheel sizes, the Start M2 Effect makes a straightforward impression. Its outstanding features include a compact folding size and durable upholstery. Thanks to its modular design the Start M2 Effect can be quickly and easily adapted to individual requirements.



Start^{M2} Effect





The patented wheel lock requires a minimum of actuating force



Individual footrests, angle-adjustable and swing-away



The rear wheel adapter can be adjusted to short or long wheelbase




The long push handles offer good grip for the attendant



The side panel can be swung completely out of the way for easy transfer



Rear wheels with quick-release axles

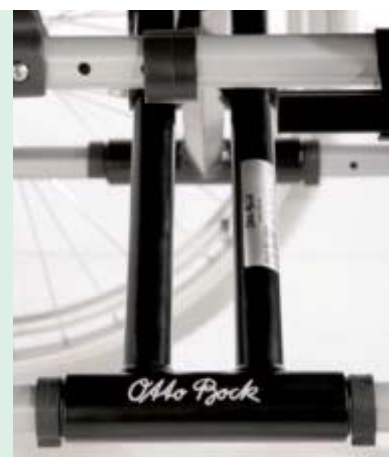

The aluminium crossbrace provides stability and good folding ability

Start^{M2} Effect



Options

Crutch holder

- Ideal for all types of walking aids



Swinging anti-tipper

- Safe in any situation



Adjustable armrest

- Also for trays
- Different versions



Start^{M2} Effect

Product Features

- Profiled aluminium tube frame
- All main component groups are attached with screws
- Standardized caster forks for all wheel sizes
- Durable upholstery
- Small folding size with folding fixture
- Stable 1 1/2-fold crossbrace
- Optimized wheel lock
- High degree of factory presetting



Technical Data

Seat width	38 cm - 48 cm	15 in - 18.9 in
Seat depth	40 cm - 46 cm	15.7 in - 18.1 in
Back height	40 cm - 47.5 cm	15.7 in - 18.7 in
Front seat height	42 cm - 53.5 cm	16.5 in - 21.1 in
Rear seat height	38 cm - 49.5 cm	15 in - 19.5 in
Lower leg length	31 cm - 47 cm	12.2 in - 18.5 in
Weight	from approx. 15 kg	from approx. 33.1 lbs
Max. user weight	125 kg	275 lbs

The Start M3 Hemi is the ideal companion for hemiplegics. Thanks to a very low anterior seat height and a wider foot space, the wheelchair user can propel the chair with his or her sound foot, provided that the footrest has been removed beforehand. Positioning shells for the affected arm as well as wheel locks for one-handed operation and one-arm drive are also available with the Start M3 Hemi.



Start^{M3} Hemi





The casters are offset outward to provide more space in the foot area



The aluminium crossbrace provides stability and good folding ability



Side frame with multiple adjustment options



Rotation adapter for modular channel forearm supports



The wheel lock can also be ordered with one-handed operation



One-arm drive with telescoping axle for independent mobility

Start^{M3} Hemi



Options

Channel forearm support

- Large selection (see Wheelchair Accessories catalogue)



Height-adjustable legrest

- Adjustable calf pad
- Continuous angle adjustment
- Automatic length adjustment



Double push ring

- Straight forward with one hand



Start^{M3} Hemi

Product Features

- Profiled aluminium tube frame
- All main component groups are attached with screws
- Standardized caster forks for all wheel sizes
- Durable upholstery
- Small folding size with folding fixture
- Stable 1 1/2-fold crossbrace
- Optimized wheel lock
- Rear wheel adapter for low anterior seat height
- Wider foot space



Technical Data


Seat width	38 cm - 50.5 cm	15 in - 19.9 in
Seat depth	40 cm - 46 cm	15.7 in - 18.1 in
Back height	40 cm - 47.5 cm	15.7 in - 18.7 in
Front seat height	38.5 cm - 50 cm	15.2 in - 19.7 in
Rear seat height	38 cm - 49.5 cm	15 in - 19.5 in
Lower leg length	31 cm - 47 cm	12.2 in - 18.5 in
Weight	from approx. 15.5 kg	from approx. 34.2 lbs
Max. user weight	125 kg	275 lbs

The Start M4 XXL is the robust model among the standard adaptive wheelchairs. With its stable double cross-brace and the reinforced rear wheels, it can carry up to 160 kilograms and withstand extreme stress and strain. An optional back stabilizer rounds out the wheelchair's features without compromising its foldability. The rear wheel adapter of the Start M4 XXL allows the centre of gravity to be adjusted to meet different user needs.



Start^{M4} XXL



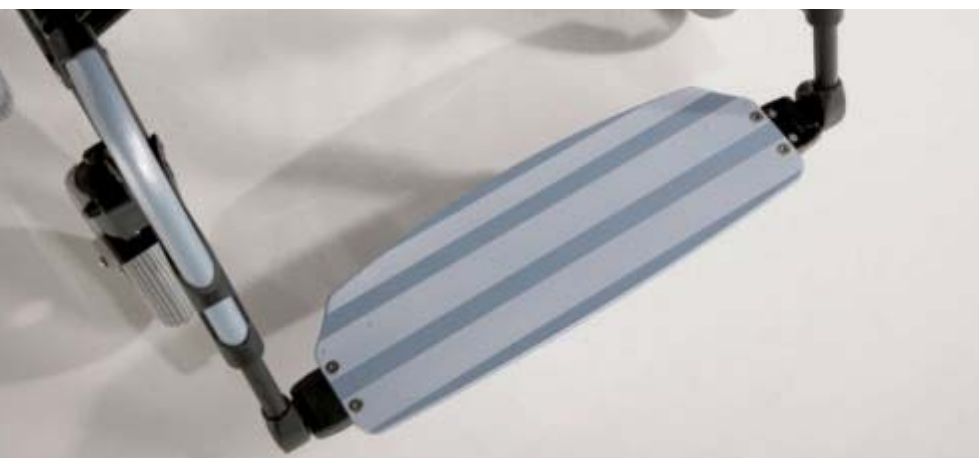

Continuous aluminium footplate for maximum space requirements



Angle adjustment of the footplate



Back stabilizer, increases overall stability and can be swung out of the way for folding



Flexible swinging anti-tipper



Rear wheel with drum brake system for attendant



Brake lever for attendant

Start^{M4} XXL



Options

Brake lever for attendant

- Also lockable



Swinging anti-tipper

- Adaptable



Adjustable armrest

- Individually adjustable



Start^{M4} XXL

Product Features

- Profiled aluminium tube frame
- All main component groups are attached with screws
- Standardized caster forks for all wheel sizes
- Durable upholstery
- Small folding size with folding fixture
- Stable double crossbrace
- Optimized wheel lock
- Reinforced rear wheels and back tubes
- Optional stabilizer bar



Technical Data

Seat width	43 cm - 58 cm	16.9 in - 22.8 in
Seat depth	43 cm - 49 cm	16.9 in - 19.3 in
Back height	40 cm - 47.5 cm	15.7 in - 18.7 in
Front seat height	40 cm - 51.5 cm	15.7 in - 20.3 in
Rear seat height	37.5 cm - 50 cm	14.8 in - 19.7 in
Lower leg length	31 cm - 47 cm	12.2 in - 18.5 in
Weight	from approx. 17.5 kg	from approx. 38.6 lbs
Max. user weight	160 kg	353 lbs

The Start M5 Comfort is constructed according to the specific wishes of the customer and offers the complete range of options. The Start M5 Comfort can be equipped with power add-on drives and can also be custom-fabricated. Its fully configurable rear wheels and precise centre of gravity settings make it similar to wheelchairs for active use. Moreover, the wheel camber setting can be changed between 0 and 2.5 degrees. Track adjustment ensures directional stability, even with increased seat inclination or wheel camber.




Start^{M5} Comfort





Rear wheel adapter for active/passive setting and wheel camber

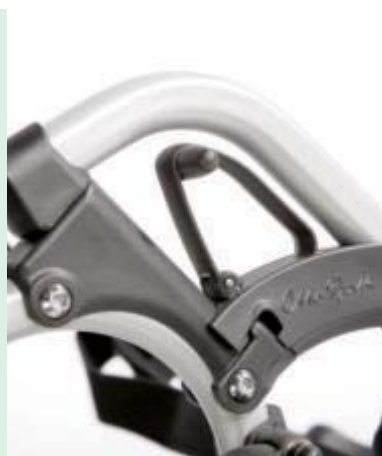

Swivel segment for footrests with quadra lock



The aluminium crossbrace provides stability and good folding ability



Height-adjustable armrest made of PU foam



Adaptable back upholstery, adjustable without the use of tools



Lightweight aluminium push ring, light-coloured anodized

Start^{M5} Comfort



Options

Plug-on clothing protector

- Individually adaptable



Amputation legrest

- Angle-adjustable



Transport wheel


- Different settings possible



Start^{M5} Comfort

Product Features

- Profiled aluminium tube frame
- All main component groups are attached with screws
- Standardized caster forks for all wheel sizes
- Durable upholstery
- Small folding size with folding fixture
- Rear wheel adapter for active/passive setting
- Optimized wheel lock
- 0° or 2.5° wheel camber
- Broad range of options

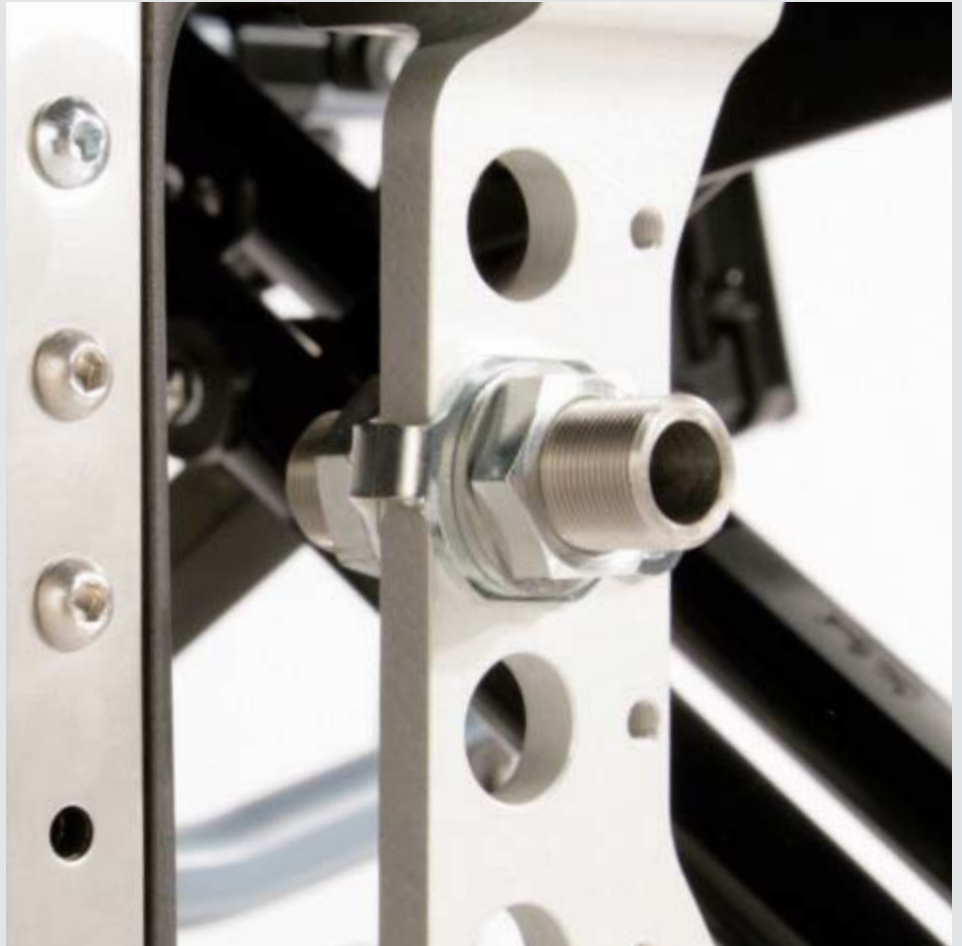

Technical Data

Seat width	38 cm - 50.5 cm	15 in - 19.9 in
Seat depth	40 cm - 46 cm	15.7 in - 18.1 in
Back height	40 cm - 47.5 cm	15.7 in - 18.7 in
Front seat height	40 cm - 51.5 cm	15.7 in - 20.3 in
Rear seat height	37.5 cm - 50 cm	14.8 in - 19.7 in
Lower leg length	31 cm - 47 cm	12.2 in - 18.5 in
Weight	from approx. 14.5 kg	from approx. 32 lbs
Max. user weight	125 kg	275 lbs

The Start M6 Junior is precisely tailored to the requirements of children. Its outstanding features include adjustable seat depth settings, clear colours, and durable upholstery. The footplate for short lower leg length as well as the casters and rear wheels were selected to match the dimensions of its users.



Start^{M6} Junior





Swivel segment for footrests with quadra lock



Standard side panel with height adjustment



The aluminium crossbrace provides stability and good folding ability



Wheel lock for user and attendant use



Mounted clamp of the anchor point kit for transportation in a motor vehicle



Rear wheels with quick-release axles

Start^{M6} Junior



Options

Height-adjustable push handle and attendant brake

- V-shape
- Swings inwards
- Removable



Rear wheel adapter allows for different settings

- From passive to active



Preparation for transportation in a motor vehicle for the disabled

- Securing points for safe anchoring in a motor vehicle for the disabled



Start^{M6} Junior

Product Features


- Profiled aluminium tube frame
- All main component groups are attached with screws
- Standardized caster forks for all wheel sizes
- Durable upholstery
- Small folding size with folding fixture
- Children's colours
- Optimized wheel lock
- Footplate for short lower leg lengths



Technical Data

Seat width	28 cm - 38 cm	11 in - 15 in
Seat depth	31 cm - 37 cm	12.2 in - 14.6 in
Back height	30 cm - 37.5 cm	11.8 in - 14.8 in
Front seat height	40 cm - 49 cm	15.7 in - 19.3 in
Rear seat height	35.5 cm - 49.5 cm	14 in - 19.5 in
Lower leg length	16 cm - 47 cm	6.3 in - 18.5 in
Weight	from approx. 14.5 kg	from approx. 32 lbs
Max. user weight	90 kg	198 lbs


Frame Colours

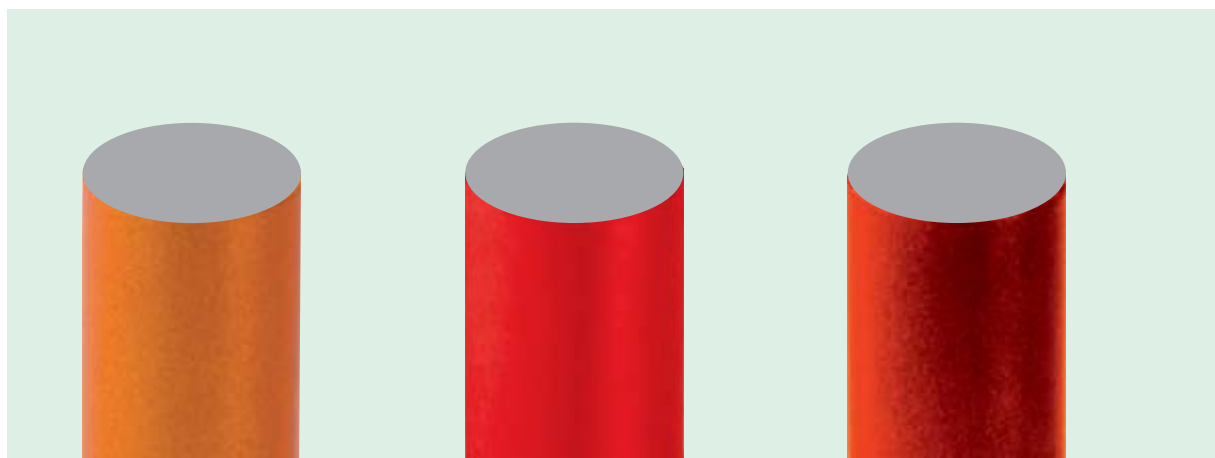

Sparkle light orange

Signal red RAL 3001

Sparkle red

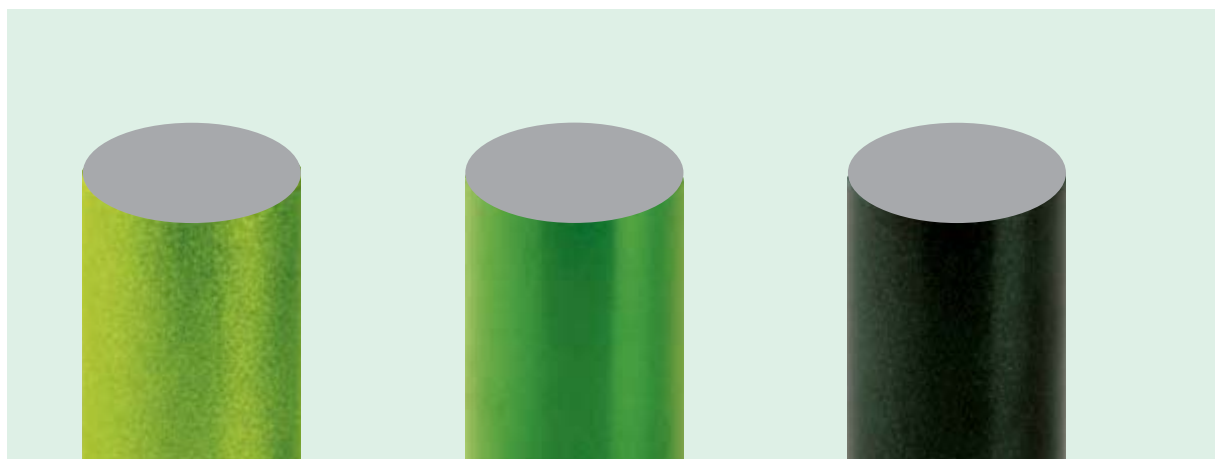

Granny smith

Dormant apple green

Sparkle blue green



Mystery

Brilliant silver

Sparkle silver new



Candy red



Dormant Vinho sparkle



Snake clear red



Clear reptile glossy



Saturn blue metallic




Dormant marine blue



Dormant sky blue



Sparkle blue

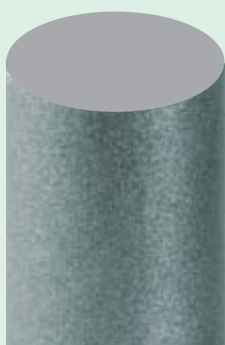

Ice-blue



Anthracite metallic



Grey matte



Deep black
RAL 9005

Notes



A large green rectangular area containing horizontal lines for writing notes. The area is filled with a light green color and contains 25 horizontal lines for writing.

[illegible]



Otto Bock®

QUALITY FOR LIFE

Otto Bock HealthCare GmbH
Max-Näder-Straße 15 · 37115 Duderstadt/Germany · Phone +49 5527 848-1590 · Fax +49 5527 848-1676 · reha-export@ottobock.de · www.ottobock.com